



QCM MATLAB Program Manual






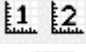

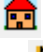



The Shull Research Group
Version 2.0b “bigfoot”

Chyi-Huey Joshua Yeh

06/17/2015

Contents

I. Quickstart Guide	3
1. Installation	3
1.1. Installing MyVNA program	3
1.1.1. Step-by-step myVNA installation	3
1.2. Installing AccessMyVNAv0.7 program	4
1.2.1. Note for code developers	5
1.3. Installing the QCM MATLAB Program	5
2. How to take measurements (Quickstart guide)	6
2.1. Starting up and closing down the programs	6
2.2. Setting up parameters and measurement options	6
2.2.1. Choosing which harmonic to track	7
2.2.2. Spectra refresh resolution	7
2.2.3. Reference time	8
2.2.4. Time between harmonics	8
2.2.5. Number of datapoints	8
2.2.6. “Record data every _ s” option	9
2.2.7. Fitting options	9
2.2.8. Formatting output file	10
2.2.9. Loading and saving settings	12
2.3. Record scan	12
2.4. Output files	12
2.4.1. No designated output filename and file location	12
2.4.2. Designated output filename and file location	14
2.5. Troubleshooting	14
II. Comprehensive guide to the QCM MATLAB Program	15
3. Features of the QCM MATLAB Program	15
3.1. Figure toolbar	15
3.1.1. Zoom in/out, pan, and data cursor tool 	16
3.1.2. Delete datapoints mode 	16

3.1.3. Load bare crystal data		16
3.1.4. Refit pre-existing raw spectras		16
3.1.5. Multi-peak fitting options		16
3.1.6. Save frequency shift data		16
3.1.7. Email notifications		16
3.1.8. Plot options		17
3.1.9. Set preferences		17
3.1.10. Home		17
3.1.11. MyVNA		17
3.1.12. Debugging mode		17
3.1.13. Help		17
3.2. GUI toolbar		18
3.3. Settings panel		18
3.4. Raw spectra panel		18
3.5. Main plot panel		18
3.6. Status toolbar		18
4. Troubleshooting		18
III. Design overview		18
5. Goal and purpose of the QCM MATLAB program		18
6. Hardware and software components		18
6.1. Hardware		18
6.2. Software		19
6.3. Design overview of how the programs work		20
IV. MATLAB code overview		22
7. Record scan algorithm explained		22
8. List of functions		22

Part I.

Quickstart Guide

1. Installation

1.1. Installing MyVNA program

The MyVNA program is a software that was developed by the original hardware designers of the N2PK Vector Network Analyzer (N2PK-VNA). This program directly manipulates the N2PK-VNA. For more information on the MyVNA program and the N2PK-VNA hardware setup, visit <http://g8kbb.co.uk/html/myvna.html>. Also, there is a small active community of radio-enthusiasts that constantly improve the software and hardware design and capabilities of the N2PK-VNA. For more information on this community, visit <https://groups.yahoo.com/neo/groups/N2PK-VNA/info>.

1.1.1. Step-by-step myVNA installation

1. **Open up the folder name, “myVNA setup files”.**
2. **Run myVNA.msi and follow the directions.** This installs the myVNA program that directly communicates and controls the N2PK Vector Network Analyzer (impedance analyzer).
3. **After installing myVNA, double click on myVNA_setup.reg.** This adds registry values associated with the myVNA program.
4. **Go to the folder in which the myVNA program is installed.** The path usually looks something like this: C:\Program Files (x86)\G8KBB\myVNA.
5. **In this folder, make sure that there is a file name, “InitGenericUsb.dll”.** If not, go to the “USB configure and Drivers” folder, located in the same path directory as the myVNA program and copy the file, “InitGenericUsb.dll” into the myVNA path directory. This file is responsible for communicating with the N2PK Vector Network Analyzer via a USB cord.
6. **In the myVNA program path directory, search for the file name, “RegServers.exe”, right click on it, and click on “Run as administrator”.** It is important that you run that execution file with administrator rights, otherwise, the execution file will not run properly. This file is responsible for registering the myVNA program with OLE capabilities. In other words, this file allows the MyVNA program to be controlled remotely with a 3rd-party program (in this case, AccessMyVNA).
7. **When myVNA program has been successfully installed, connect the N2PK VNA to the computer.** The computer will automatically start scanning for the appropriate drivers to download and install, in order to interface with the N2PK VNA. It is not necessary for the computer to go through this process (canceling this process will not affect the driver installation), since the associated driver will be installed manually. Go to Control Panel > View devices and printers. Scroll down and right click “Unknown device”. From there, the driver can be manually installed by directing Windows to look into the

folder in which the myVNA program was installed. Make sure “Check subfolder” option is checked.

8. When running myVNA for the first time, a message will display, saying that the “Calibration” file was not found. Click “OK”; the GUI will initialize. The first thing to do is to calibrate the N2PK VNA. Go to “Calibration Options” (left panel) and turn on the “Unguided Calibration” option. On the top toolbar, set the start and end frequency to 4 MHz and 56 MHz. Make sure the units are set to “MHz” and that the “Start/Stop” setting is selected. Also, change the number of setps to 400 and the average to 3. Click on the “Calibrate” button, located at the top toolbar.

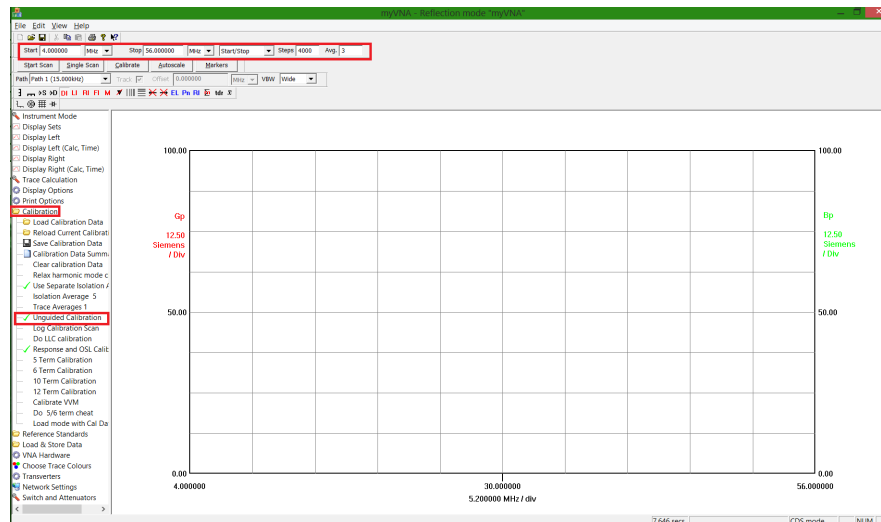


Figure 1: Set up te calibration process.

9. Make sure that the N2PK VNA is in an open-circuit state. Click on the “Open” button. The program will perform a scan from 4 MHz to 56 MHz with 4000 datapoints in between. Once this is finished, connect the “short” to the N2PK VNA, and click on the “Short” button. Repeat this step with the 50 Ω resistor and click on the “Load” button. Once the calibration is finished, click on the the “Finished” button and save the calibration file.
10. The myVNA program should fully functioning. To test it, simply click on the “Single Scan” button on the top toolbar. The “Autoscale” button might need to be pressed in order to adjust the plot to show the measurement. Once everything has been tested and properly set-up, exit out of myVNA.

Note: It is good practice re-calibrate before an experiment. This will reduce any unnecessary complication associated with the Lorentzian fitting process QCM measurement.

1.2. Installing AccessMyVNAv0.7 program

Setting up this program is easy. All that is required is to ensure that the folder, “AccessMyVNAv0.7”, is in the same directory as the QCM MATLAB Program. For example, if the path directory to the QCM MATLAB Program is <C:\Program Files (x86)\QCM MATLAB Program>, make sure that <C:\Program Files (x86)\QCM MATLAB Program\AccessMyVNAv0.7> exists. This folder contains the program, AccessMyVNA, that is important for the QCM MATLAB Program to (indirectly) manipulate the MyVNA program. Currently, the QCM MATLAB Program automatically executes the AccessMyVNA each time the QCM MATLAB Program initiates

and creates the GUI figure. However, to manually execute the AccessMyVNA program, the “AccessMyVNA.exe” file can be found in the “release” folder in the AccessMyVNAv0.7 folder.

Make sure to install the necessary .dll files. To do this, simply run execution file located in the “Visual C++ redistributable packages for visual basic 2013” folder. If the installation is being performed on a 32-bit computer, run the “vcredist_x86(1).exe” file. If the installation is being performed on a 64-bit computer, run the “vcredist_x64(1).exe” file. If an error message appears, (when executing the AccessMyVNA.exe file) stating that a .dll file is missing, try running both the 32-bit and 64-bit execution files.

For those running the “AccessMyVNAv0.7” on Windows 8 or higher, make sure to change the compatibility mode of the “AccessMyVNAv0.7.exe” to Windows 7. Otherwise, the program will unexpectedly crash.

1.2.1. Note for code developers

Within the AccessMyVNAv0.7 folder, there are many other folders and files. These “extra” folders and files are needed if there is a need to edit the AccessMyVNAv0.7 program code. To access this code, Visual Basic C++ (2013 or higher) must be installed. The file, “AccessMyVNA.sln”, can be opened with Visual Basic C++ and the file, AccessMyVNADlg.cpp, can be accessed in the solution explorer panel. It is within the AccessMyVNADlg.cpp in which most of the important features and functions of the AccessMyVNAv0.7 is defined. Ideally, this program should be rewritten. However, rewriting this program requires advanced knowledge and experience in working with Visual Basic C++. In particular, this program relies on “MFC” (Microsoft Foundation Class) libraries, which are only supported in Visual Basic Professional 2013 or higher (not the “Express” version!). A free student license of Visual Basic Professional 2013 (or higher) can be obtained from Microsoft’s Dreamspark program.

1.3. Installing the QCM MATLAB Program

The QCM MATLAB Program does not require any formal installation. However, make sure the program is performed by MATLAB 2014b or higher. Make sure that the following files:

1. “QCM_v002b_bigfoot.m”
2. “QCM_v002b_bigfoot.fig”
3. “fg_values.mat”
4. “raw_spectras.mat”

are in the <C:\Program Files (x86)\QCM MATLAB Program\> path directory. Also, make sure the following folders are located in the <C:\Program Files (x86)\QCM MATLAB Program\> path directory:

1. “QCM MATLAB manual”
2. “email files”
3. “AccessMyVNAv0.7”
4. “deleted data”

Note that the filenames and folder names are case sensitive. To run the QCM MATLAB Program, double click on the “QCM_v001a.m” file and run the script. As aforementioned, the AccessMyVNA program should also start with the creation of the MATLAB GUI figure. If it did not, an error should be reflected in the MATLAB command window.

2. How to take measurements (Quickstart guide)

2.1. Starting up and closing down the programs

To start up the QCM MATLAB program (and the other programs, a.k.a. MyVNA and AccessMyVNA), run the “QCM_v001a.m” script. An instance of the AccessMyVNA program should appear along with the creation of the MATLAB GUI figure. In the AccessMyVNA window, click on the button, “Start Scan” (see Figure 2). This will automatically start up the MyVNA program and begin (continuous) scans. Note that when the AccessMyVNA program is continuously scanning, the program cannot be minimized and may seem unresponsive. As long the AccessMyVNA program is telling the myVNA program to continue scanning, the AccessMyVNA program is still running correctly.

To have AccessMyVNA program to stop scanning, simply go to the radio dial labelled, “Maintain myVNA scan”, and uncheck the radio dial. This will communicate to the AccessMyVNA program to stop scanning (this process may take a few seconds to complete). If the program does not stop scanning, then the AccessMyVNA program is probably malfunctioning. In that case, the AccessMyVNA program may have to be terminated in the windows task manager. **To begin the continous scanning again, make sure that the radiodial, “Maintain myVNA scan” is checked and click on the “Start Scan” button in the AccessMyVNA window.** Also, note that when the MATLAB GUI figure is closed, the QCM MATLAB code will communicate the AccessMyVNA to stop scanning. However, the **AccessMyVNA program and the MyVNA program needs to be manually closed.** If the AccessMyVNA program is not manually closed, reinitiating the QCM MATLAB GUI figure will create a new instance of AccessMyVNA; thus, there will be multiple instances of AccessMyVNA running at the same time, which can cause fatal errors.

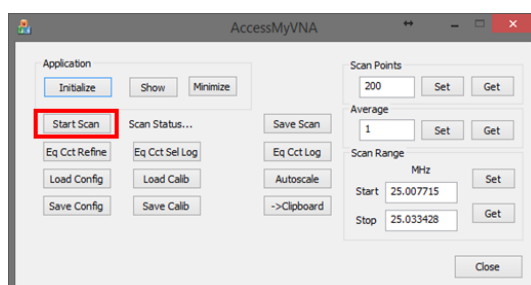


Figure 2: Click on the “Start Scan” button in the AccessMyVNA program to begin scanning.

2.2. Setting up parameters and measurement options

Parameters, settings, and options related to the collecting scans from the N2PK-VNA impedance analyzer can be controlled from the settings panel (see Figure 3).

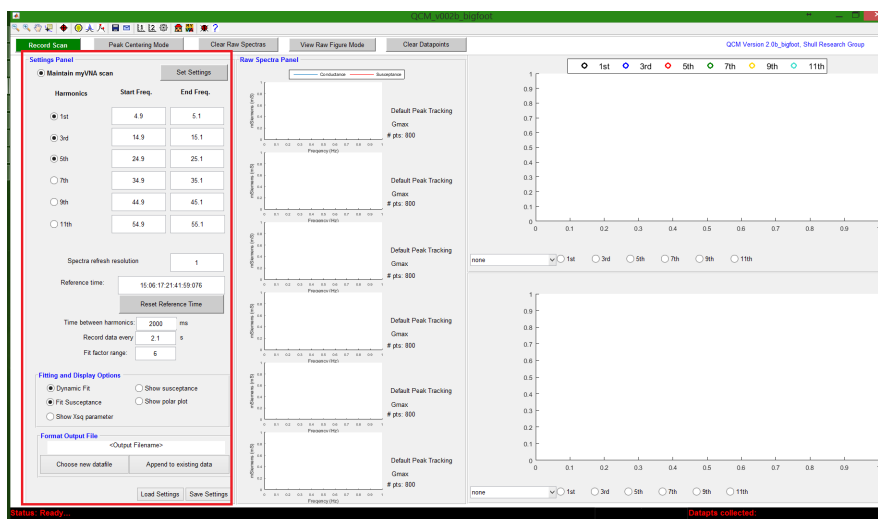


Figure 3: The settings panel contains information on the parameters, settings, and options related to the collection of the raw data.

2.2.1. Choosing which harmonic to track

The N2PK-VNA impedance analyzer is capable of tracking up to 6 harmonics at the same time (1st, 3rd, 5th, 7th, 9th, 11th harmonics). To choose which harmonics to track, check/uncheck the radio dials on the “Harmonics” column (see Figure). The start and end frequencies of the scan for each of the harmonic can be set manually by editing the numbers in the corresponding start/end frequencies. Also, the start and end frequencies can be set by initiating the “Peak Centering Mode”. For more details in regards to the Peak Centering Mode, refer to Part 1, Section 3.1.

Harmonics	Start Freq.	End Freq.
<input checked="" type="radio"/> 1st	4.9	5.1
<input type="radio"/> 3rd	14.8	15.1
<input type="radio"/> 5th	24.8	25.1
<input type="radio"/> 7th	34.9	35.1
<input type="radio"/> 9th	44.9	45.1
<input type="radio"/> 11th	54.9	55.1

Figure 4: Multiple harmonics can be tracked. Choose which harmonics to track by checking/unchecking the radio dials under the “Harmonics” column. The start and end frequencies of the scan can be set for each harmonic.

2.2.2. Spectra refresh resolution

This option in the “Settings Panel” controls how often the QCM MATLAB Program refreshes the raw conductance and/or susceptance spectras. Its default value is 1, which means that the MATLAB program will refresh the spectras each time it reads in data from the datafile (see Part 2 for a more thorough explanation on how the QCM MATLAB Program reads data from the impedance analyzer). Changing this number to a higher integer-value can potentially increase the time resolution of the measurements significantly, since it takes a finite amount of time for MATLAB to constantly refresh the spectral plots. For example, a spectra refresh resolution of 10 means that the MATLAB program will refresh the plots every 10th measurement it collects.

2.2.3. Reference time

Since this programs collects frequency and dissipation shifts as a function of time, having a reference time is important, especially if new data needs to be appended onto another dataset from a previous experiment. The default reference time is set to the time in which the QCM MATLAB GUI figure was initiated. The reference time can also be changed manually by editing the values in the text box. Make sure that the values are in the correct format, yy:mm:dd:HH:MM:SS:FFF, where yy, mm, dd, MM, SS, and FFF, represent the last two digits of the year, the month, the day, the hour (24-hour format), the minutes, the seconds, and the milliseconds, respectively. **Do not hit the “Reset Reference Time” button after editing the reference time. Clicking on the “Reset Reference Time” will set the reference time to the current time! Simply hit the “Enter” key after inputting your desired reference time.**

2.2.4. Time between harmonics

This option controls the amount of time (in milliseconds) in between each harmonic measurements. In general, it is not necessary to change this value. **A short time can lead to synchronization and systematic errors.** If the resonance peaks in the spectra looks odd or abnormal, increasing the amount of time in between each harmonic measurement might solve the problem. Otherwise, it might be due to other errors (see the Troubleshooting section). An example of what a resonance peak looks like if there is not sufficient amount of time in between harmonics is shown in Figure 5.

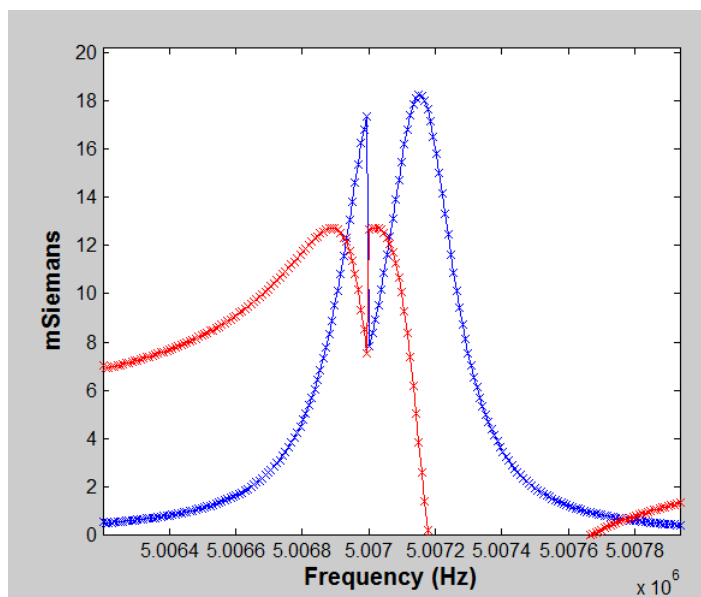


Figure 5: An abnormal resonance peak (5MHz) due to insufficient amount of time in between harmonics. Increase the time in between each harmonic measurement to solve this problem.

2.2.5. Number of datapoints

The number of datapoints in each frequency scan can be adjusted. The default value is 200 datapoints. Note that increasing the number of datapoints will require more time in between harmonics. Thus, the time in between harmonics need to be adjusted accordingly. **A good rule of thumb is to allow 4-6 ms for each datapoint.**

2.2.6. “Record data every _ s” option

The time interval in which the QCM MATLAB Program records the spectra measure can be set in this option. For experiments that last for days to weeks, it is more practical to increase the time interval in which the MATLAB program records the data. Currently, the MATLAB program allows for 1 million timepoints to be recorded. If more than one million datapoints are collected, an error will appear and terminate the data collection process. The maximum number of datapoints to keep can be adjusted in the MATLAB code. However, it should be noted that the MATLAB program will progressively slow down as more data is being collected since it will take longer for MATLAB to append new values and variables into the output “.mat” files. For all practical purposes, this should not be an issue. If it is an issue, the output files should be split up into multiple files to mitigate this effect.

For experiments that are on a short time scale, it is important to remember that the number of datapoints and time in between harmonics need to be adjusted accordingly before decreasing the recording time interval. Otherwise, duplicates of the same spectra will be recorded. **In other words, if the recording time interval is less than the time in between harmonics, duplicates of the same spectra will be recorded.** An example of this systematic error is shown in Figure 6.

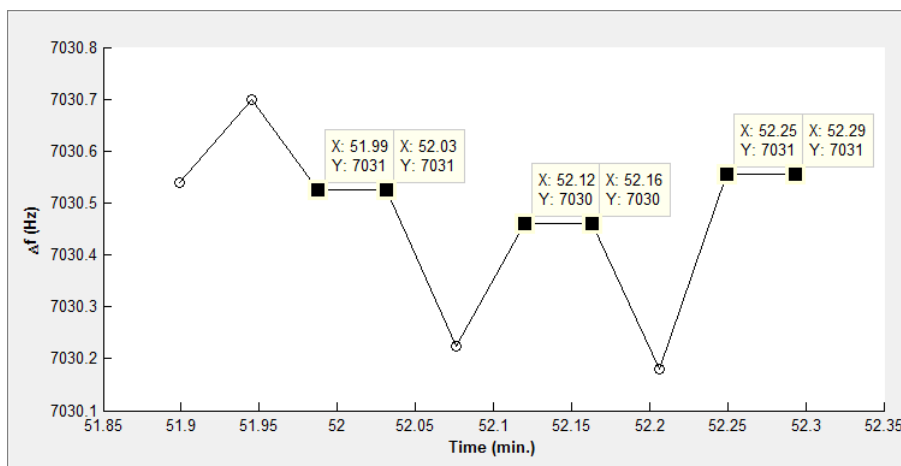


Figure 6: If the recording time interval is less than the time in between harmonics, duplicates of datapoints will occur. This is evident if the frequency shift plot is observed during the measurement process.

2.2.7. Fitting options

There are five radio-dials in the “Fitting Options” panel (see Figure 7): “Dynamic Fit”, “Fit Susceptance”, “Show Xsq parameter”, “Show susceptance”, and “Show polar plot”. If the “Dynamic Fit” radiodial is checked, a Lorentzian function curve fitting algorithm will be used to fit the resonance conductance peaks each time a spectra is collected. Unchecking this option significantly decreases the time required to record a measured spectra, which may be of importance if recording fast scans are required in an experiment. If the “Fit Susceptance” radiodial is checked along with the “Dynamic Fit” radio dial, both the conductance and susceptance resonance curves will be fitted. Note that the additional step of fitting the susceptance curve will require more time to record a measured spectra. If the “Show Xsq parameter” radiodial is checked along with the “Dynamic Fit” radiodial, the χ^2 parameter for each fit will be calculated and recorded. Details in regards to calculating goodness-of-fit for the curve fitting process is discussed in later sections.

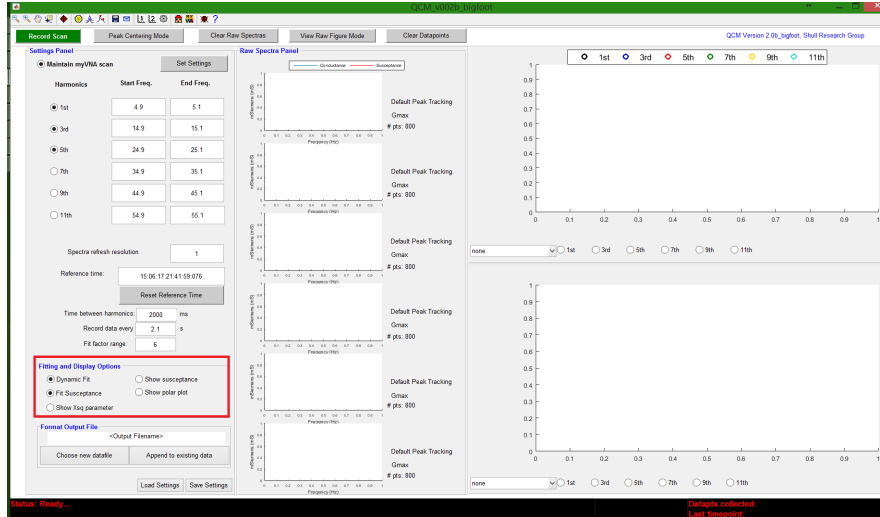


Figure 7: Options related to the fitting process during the measurements can be adjusted in the “Fitting Options” section located in the “Settings” panel.

2.2.8. Formatting output file

In the “Format Output File” panel, there is a button called, “Choose data filename”. This button will open a file-output explorer dialogue box (see Figure 8). From this dialogue box, the name and location of the output data can be chosen. **If a pre-existing file is chosen, any absolute frequency/dissipation and frequency/dissipation shift data that is collected will be overwritten unless the settings associated with the output datafile is also loaded.** If the frequency/dissipation data was not recorded (in other words, a Lorentzian function was not fitted to the raw spectra data), the saved raw spectra curves will be added into the .mat file containing the raw spectras. **However, be sure to set the correct reference time associated with the original output datafile. Otherwise, the timestamp associated with the raw spectras will be incorrect.** Generally speaking, choosing a pre-existing file is not recommended. Appending new data to a pre-existing file can be better accomplished by clicking on the “Append existing data” button. This allows the user to load pre-existing data sets. Data collection will added to the pre-existing datafiles. For more information on how the data is stored in the output data file, see the section on “Loading and saving settings”.


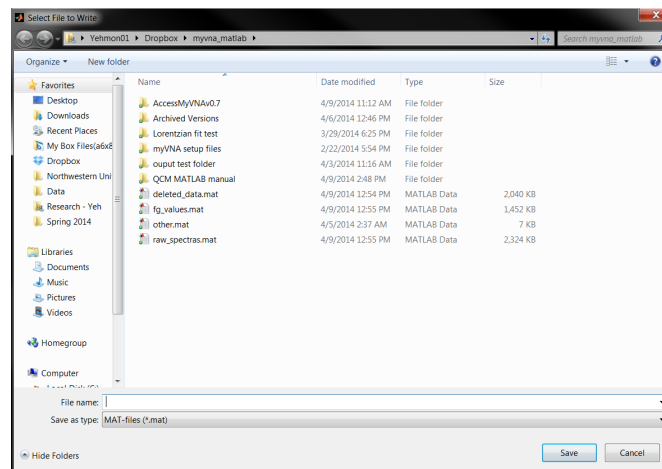



Figure 8: File-output explorer dialogue box.

If the “settings” file associated with the output datafile does not exist, the frequency/dissipation data can still be manually appended. To do this, follow the steps listed below:


1. Click on the “Access handles structure” button.

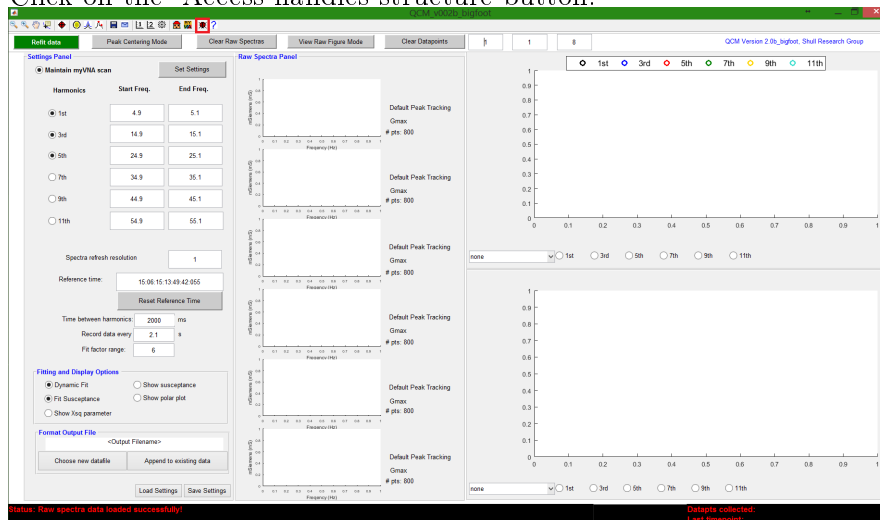

2. Go to the MATLAB Command Window and push Ctrl+o buttons on the keyboard.

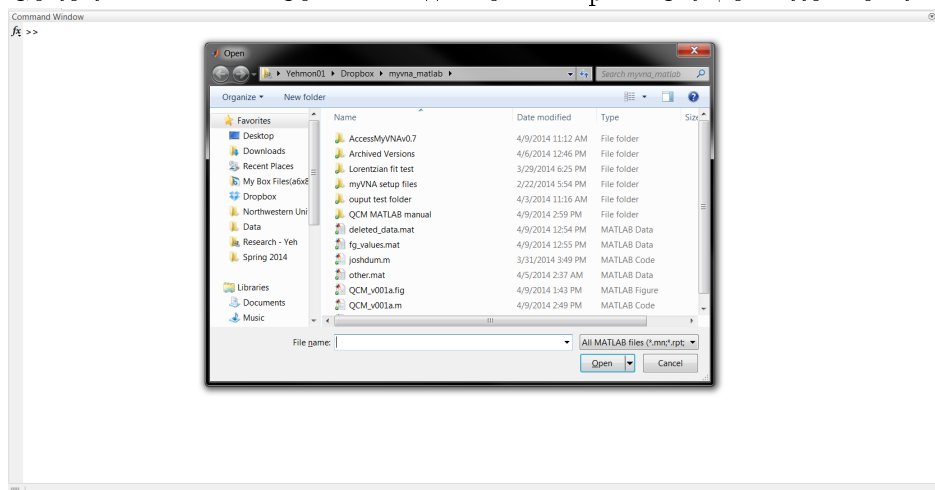


3. Open the output datafile that was chosen previously from the file-output explorer dialogue. This will load the variables stored in the .mat file into the MATLAB Workspace.
4. In the MATLAB Workspace, double-click on abs_freq or freq_shift

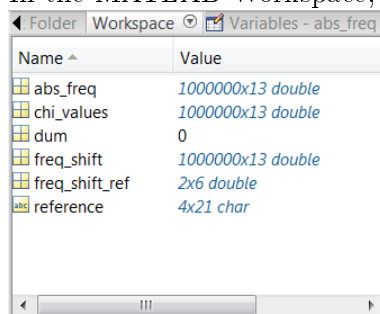

5. This will open up the “Variables” dialogue box containing the data stored in the variable that was double-clicked. Find the row number associated with the last datapoint that was stored.

- Go back to the Command Window and type in the following command:
“handles.din.n=<row number associated with the last datapoint that was stored>”.
This will tell Matlab where (based on the row number) to begin appending the new frequency and dissipation data.
- Click into the Matlab Editor Window and click on the “Quit Debugging” button near the top of the window.

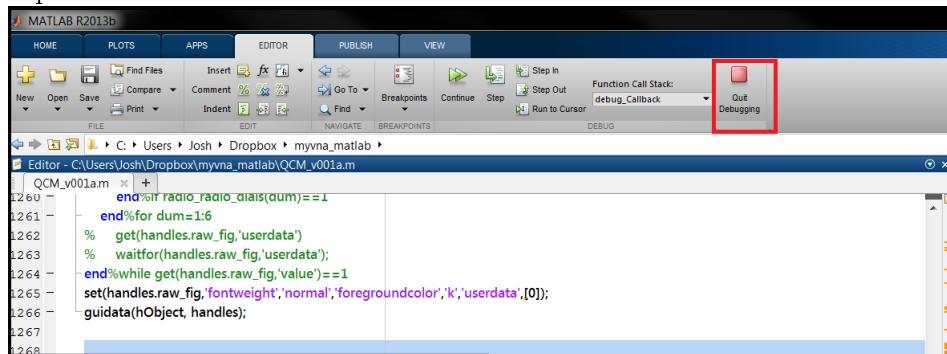


2.2.9. Loading and saving settings

The options and setting parameters can be saved. By clicking the “Save Settings” button, a settings file will be saved in the same path location as the output datafiles. If no output datafile was selected to begin with, the the current setting options and parameters will be saved in the same path directory as the QCM MATLAB Program.

Clicking on the “Load Settings” button will automatically load the settings file associated with output datafile (which is programmed to be named as “<Output Filename>_settings.mat”). If the automatically-loaded settings file was not the desired settings to be loaded, the “Load Settings” button can be clicked again to manually load the desired settings file. If the settings file cannot be found, an open file dialogue box will pop-up and the desired settings file can be loaded.

When the QCM MATLAB Program runs into any errors while loading or saving the setting options and parameters, the program will go into “Debugging mode”. To exit out of the debugging mode, simply push the button that says, “Quit debugging”, or type “return” in the MATLAB Command Window.

2.3. Record scan

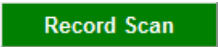To begin recording the measurements, simply click on the “Record scan” button . Details on how the measurement is recorded can be seen in the QCM MATLAB code. That section of the code is thoroughly commented. Also, refer to later sections on the summary of how the MATLAB program records and saves the measurement scans.

2.4. Output files

2.4.1. No designated output filename and file location

When no output file is designated in the “Format Output File” panel in the “Settings” section, the recorded data will be stored in the default output files:

- “fg_values.mat”: This MATLAB .mat file contains frequency and dissipation information that was collected during the recording process.

Variables

- a) **“abs_freq”**: (1e6)x13 double array containing the absolute frequency and dissipation data. The first column contains the timepoints. Columns 2, 4, 6, 8, 10, and 12 contain the absolute frequency data. Columns 3, 5, 7, 9, 11, and 13 contain the absolute dissipation data.
 - b) **“freq_shift”**: (1e6)x13 double array containing the frequency and dissipation shift data. The first column contains the timepoints. Columns 2, 4, 6, 8, 10, and 12 contain the frequency shift data. Columns 3, 5, 7, 9, 11, and 13 contain the dissipation shift data.
 - c) **“chi_values”**: (1e6)x3 double array containing the χ^2 curve fitting statistics. The first column contains the timepoints. Columns 2 and 3 contain the χ^2 curve fitting statistics for the conductance and susceptance data.
 - d) **“freq_shift_ref”**: contains the reference frequency (row 1) and dissipation (row 2) values for the frequency and dissipation shift data. Columns 1, 2, 3, 4, 5, and 6 contain reference values for harmonics 1, 3, 5, 7, 9, and 11.
 - e) **“reference”**: 4x1 cell array containing (in order of row index) the reference timestamp (yy:mm:dd:HH:MM:SS:FFF), the time in between harmonic measurements in ms, the number of datapoints collected in the raw spectra data, and the time interval in which the data was recorded in s.
 - f) **“version”**: string variable containing information on which version of the QCM program was used to record the data
2. **“raw_spectras.mat”**: This MATLAB .mat file contains raw spectra information that was collected during the scan.

Variables

- a) **“reference”**: 4x1 cell array containing (in order of row index) the reference timestamp (yy:mm:dd:HH:MM:SS:FFF), the time in between harmonic measurements in ms, the number of datapoints collected in the raw spectra data, and the time interval in which the data was recorded in s.
- b) **filename_t_min.dotfractional min_iq_1_ih_harmonic_order**: contains the frequency (first column), conductance (second column), and susceptance data (third column) at time <min>.<fractional min>.

Both files are located in the same path directory as the QCM MATLAB Program. It is important to note that these two .mat files may or may not be cleared before the data collection process. Thus, any data stores may or may not be convoluted from previous stored data in the .mat files. These two files is just a safeguard from losing data that was collected. If data is extracted from these two files, it is important to clear the data that was stored, so that future data stored in the files will not be convoluted with old data. To do this, follow these steps:

1. Go to the MATLAB Command Window and type, “clear all”. Hit the “Enter” key.
2. In the MATLAB Command Window, type “dum=0”. Hit the “Enter” key.
3. Click into the MATLAB Workspace and push “Ctrl+o” buttons on the keyboard.
4. This will create a save prompt window. Save and replace both “fg_values.mat” and “raw_spectras.mat”.

If the settings are saved with no designated output file, a “default_settings.mat” file will be saved in the same path directory as the QCM MATLAB Program.

2.4.2. Designated output filename and file location

When an output filename and location are designated, the recorded data will be stored in the following output files:

1. “<user-designated filename>.mat”: This MATLAB .mat file contains frequency and dissipation information that was collected during the recording process.


Variables

- a) “abs_freq”: (1e6)x13 double array containing the absolute frequency and dissipation data. The first column contains the timepoints. Columns 2, 4, 6, 8, 10, and 12 contain the absolute frequency data. Columns 3, 5, 7, 9, 11, and 13 contain the absolute dissipation data.
 - i. “freq_shift”: (1e6)x13 double array containing the frequency and dissipation shift data. The first column contains the timepoints. Columns 2, 4, 6, 8, 10, and 12 contain the frequency shift data. Columns 3, 5, 7, 9, 11, and 13 contain the dissipation shift data.
 - ii. “chi_values”: (1e6)x3 double array containing the χ^2 curve fitting statistics. The first column contains the timepoints. Columns 2 and 3 contain the χ^2 curve fitting statistics for the conductance and susceptance data.
 - iii. “freq_shift_ref”: contains the reference frequency (row 1) and dissipation (row 2) values for the frequency and dissipation shift data. Columns 1, 2, 3, 4, 5, and 6 contain reference values for harmonics 1, 3, 5, 7, 9, and 11.
 - iv. “reference”: 4x1 cell array containing (in order of row index) the reference timestamp (yy:mm:dd:HH:MM:SS:FFF), the time in between harmonic measurements in ms, the number of datapoints collected in the raw spectra data, and the time interval in which the data was recorded in s.
2. “<user_designated filename>_raw_spectras.mat”: This MATLAB .mat file contains raw spectra information that was collected during the scan.

Variables

- a) “reference”: 4x1 cell array containing (in order of row index) the reference timestamp (yy:mm:dd:HH:MM:SS:FFF), the time in between harmonic measurements in ms, the number of datapoints collected in the raw spectra data, and the time interval in which the data was recorded in s.
 - b) filename_t_<min>dot<fractional min>_iq_1_ih_<harmonic order>: contains the frequency (first column), conductance (second column), and susceptance data (third column) at time <min>.<fractional min>.
3. If a settings file was saved, “<user-designated filename>_settings.mat”

2.5. Troubleshooting

The QCM MATLAB Program code is thoroughly commented and tested. If errors occur, the code can be examined and debugged. The easiest way to access the QCM MATLAB GUI handles structure is by pushing the  button.

Part II.

Comprehensive guide to the QCM MATLAB Program

3. Features of the QCM MATLAB Program

This section provides a more comprehensive overview of all the features available in the MATLAB program. The basic sections of the GUI are shown in .

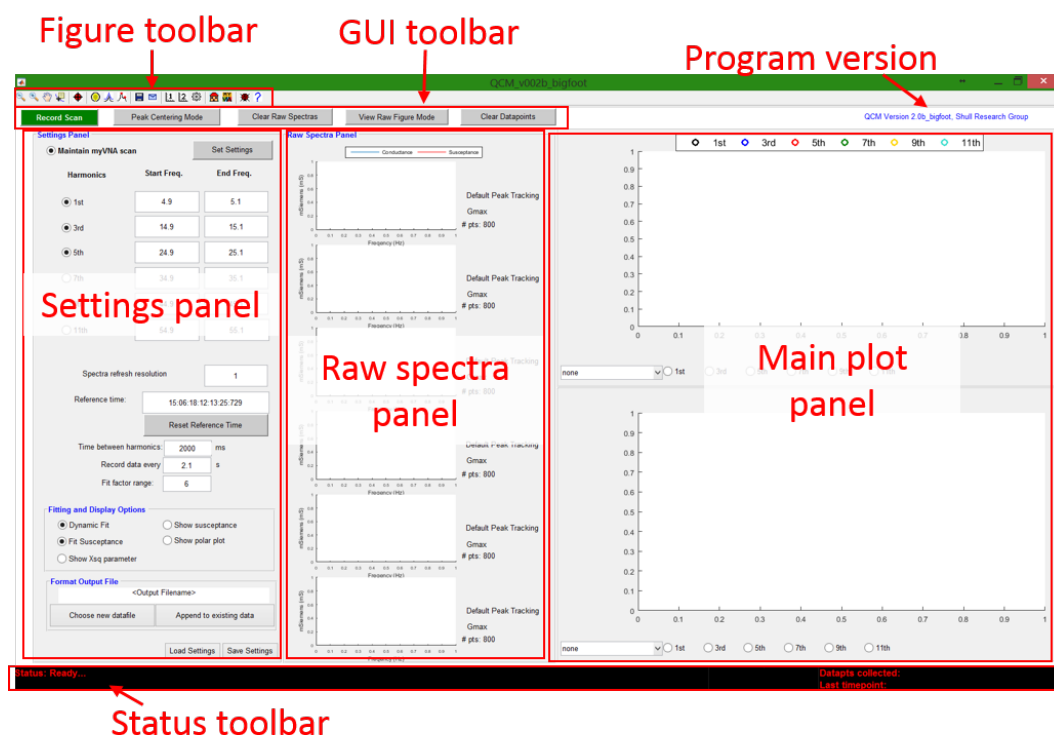


Figure 9: GUI Overview

3.1. Figure toolbar

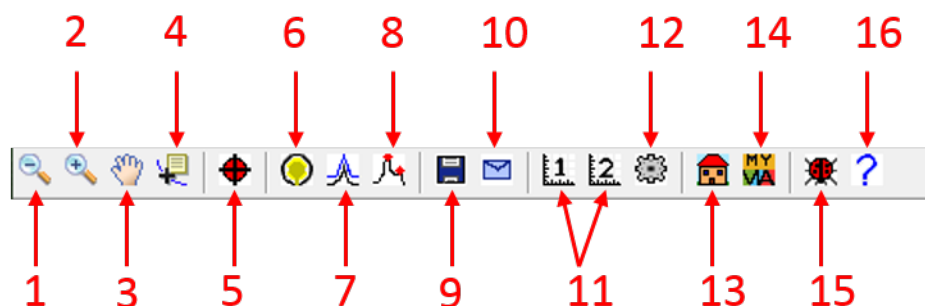


Figure 10: A list of all of the buttons on the figure toolbar

#	name	description
1	zoom-in	zoom-in into an axes plot
2	zoom-out	zoom-out from an axes plot
3	pan tool	pan axes tool
4	data cursor	select datapoints in an axes
5	delete datapoints mode	enable “delete datapoints mode”
6	load bare crystal data	load reference frequencies
7	refit pre-existing raw spectras	this will enable the “refit mode”
8	multi-peak fitting options	set multi-peak functions
9	save frequency shift data	save the frequency shift data
10	email notifications	set-up email notifications
11	plot options	plot options for primaryaxes1 and primaryaxes2
12	set preferences	set-up GUI preferences
13	home	revert GUI state back to the default state
14	MyVNA	execute AccessMyVNA.exe
15	debugging mode	enter into debugging mode (for advanced users only!)
16	help	opens the pdf manual

Table 1: Short description of each button on the figure toolbar

3.1.1. Zoom in/out, pan, and data cursor tool

These tools have the same functionality as the zoom in/out, pan, and data cursor tools on a standard figure window.

3.1.2. Delete datapoints mode

3.1.3. Load bare crystal data

3.1.4. Refit pre-existing raw spectras

3.1.5. Multi-peak fitting options

3.1.6. Save frequency shift data

3.1.7. Email notifications

Clicking on this button will bring up a window asking the user to input an email address and an option to turn on/off the email notifications.

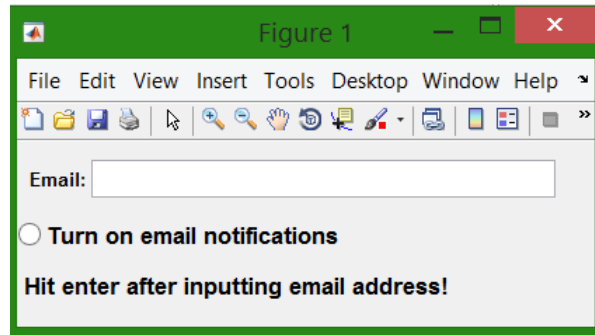


Figure 11: Email notification window.

This function only stores the email information and the toggle state of the email notifications in the handle structure. When the notifications are turned on, an email will be sent (using the “email_send” function) during the measurement process. The “email_send” function utilizes undocumented Javascript code. A screenshot of the GUI (.jpg) will be sent along with a Command Window log text file, QCM_diary.txt (see “Diary” function in the MATLAB help document). Note that the QCM_diary.txt file will continue to grow larger as more points are being recorded (since there are more “stuff” in the command window). Thus, it is advisable to clear out the QCM_diary.txt file ever once in awhile during a long measurement process with short measurement intervals. Also, if the “email_send” functions runs into an error, an error message will appear in the command window, but the measurement will continue (the main code is written in a try/catch-block).

3.1.8. Plot options



This button allows the user to adjust the properties for “primaryaxes1” (top right axes) and for “primaryaxes2” (bottom right axes). The callback function uses the “inspect” function to provide the user access to the axes properties. For more information on what each axes property does, refer to the MATLAB help document on the “inspect” function.

3.1.9. Set preferences



3.1.10. Home



3.1.11. MyVNA



3.1.12. Debugging mode



This takes the user into debugging mode. The callback function uses the “keyboard” function and also exports the handle structure of the GUI into the “base” workspace. Note that any changes that are made will not necessarily be saved. For example, changes to variables (not handle objects) will not be saved in the handles structure unless the “guidata” function has been initiated in the command window. However, direct changes to handle objects (using the “set” function), will be immediately applied and saved. For more information on the details on how information is saved, see the “guidata” in the MATLAB general help files.

3.1.13. Help



This button will direct the user to the pdf manual. Also, feel free to contact me through email (CHJoshuaYeh@u.northwestern.edu).

3.2. GUI toolbar

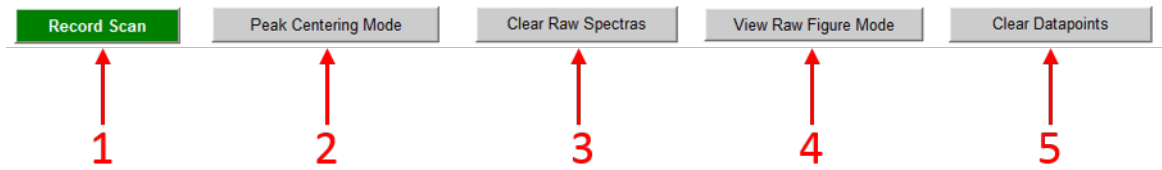


Figure 12: GUI toolbar

#	name	description
1	Record Scan	Start the measurement process
2	Peak Centering Mode	Allows for finding the appropriate harmonics
3	Clear Raw Spectras	Clear all of the spectral plots
4	View Raw Figure Mode	View each spectra in individual plots
5	Clear Datapoints	Delete all frequency shift datapoints

Table 2: Short description of each button on the GUI toolbar

3.3. Settings panel

3.4. Raw spectra panel

3.5. Main plot panel

3.6. Status toolbar

4. Troubleshooting

Part III.

Design overview

5. Goal and purpose of the QCM MATLAB program

The primary goal for the development of this program is to provide a MATLAB interface in obtaining QCM measurements. Specifically, this program has the capability to measure the location and half-max-half-width (HMH) of the resonance conductance peaks. Based on this information, analysis can be performed (in MATLAB) on these measurements to extract out the viscoelastic properties of a film deposited on the QCM crystal. Currently, error bars are not explicitly calculated (only χ^2 values for the peak fitting is calculated); however, one can easily calculate error statistics since all of the raw conductance/susceptance spectra data is saved in a “.mat” .

6. Hardware and software components

6.1. Hardware

The hardware required to run the program is the N2PK Vector Network Analyzer (N2PK-VNA, an impedance analyzer custom built by Ivan Markarov), a crystal holder (Inficon, NY), and a


QCM crystal (Inficon, NY). Further information on the N2PK-VNA hardware details can be found in <http://g8kbb.co.uk/html/myvna.html>. An image of the setup can be seen in Figure 13.

Figure 13: An image of the overall hardware setup.


6.2. Software

The N2PK-VNA came with a proprietary software called myVNA. This software allows for remote access, which provides a means to control the N2PK-VNA hardware with a custom program. Currently, the program called AccessMyVNA (written in VB C++) is used to remote-access the impedance analyzer. AccessMyVNA is designed as a “middleman” or a “gateway” program for the QCM MATLAB program to communicate with the N2PK-VNA hardware. Details on how the programs communicate with each other is discussed in subsequent sections. To summarize, in order to collect QCM measurements with the QCM MATLAB program, three programs must be running simultaneously: 1) MyVNA, 2) AccessMyVNA, 3) QCM MATLAB program. An image if each program is shown in Figure.

a) MyVNA

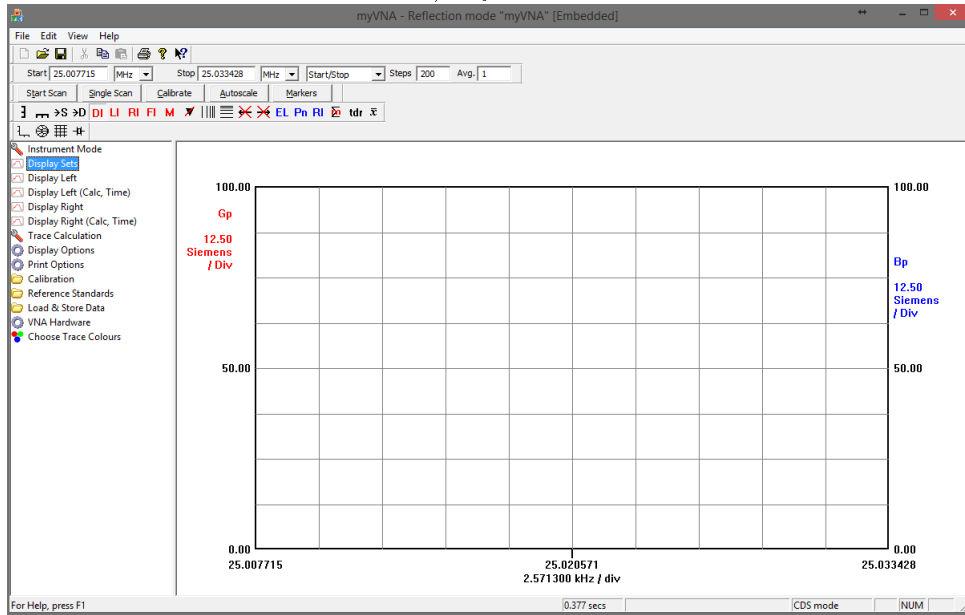

b) AccessMyVNA

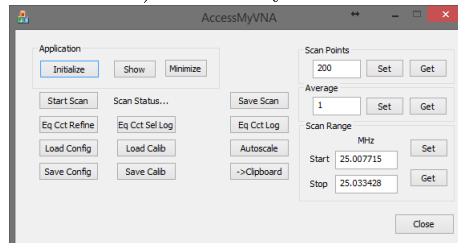

c) QCM MATLAB Program

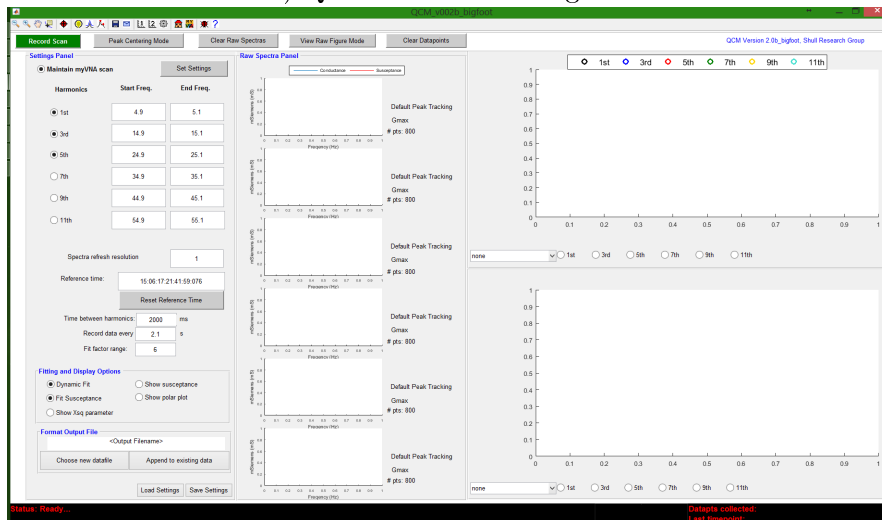


Figure 14: A screen shot of each program: a) MyVNA, b) AccessMyVNA, c) QCM MATLAB Program

6.3. Design overview of how the programs work

A general overview of how all the programs communicate with each is shown in flow chart in Figure 15.

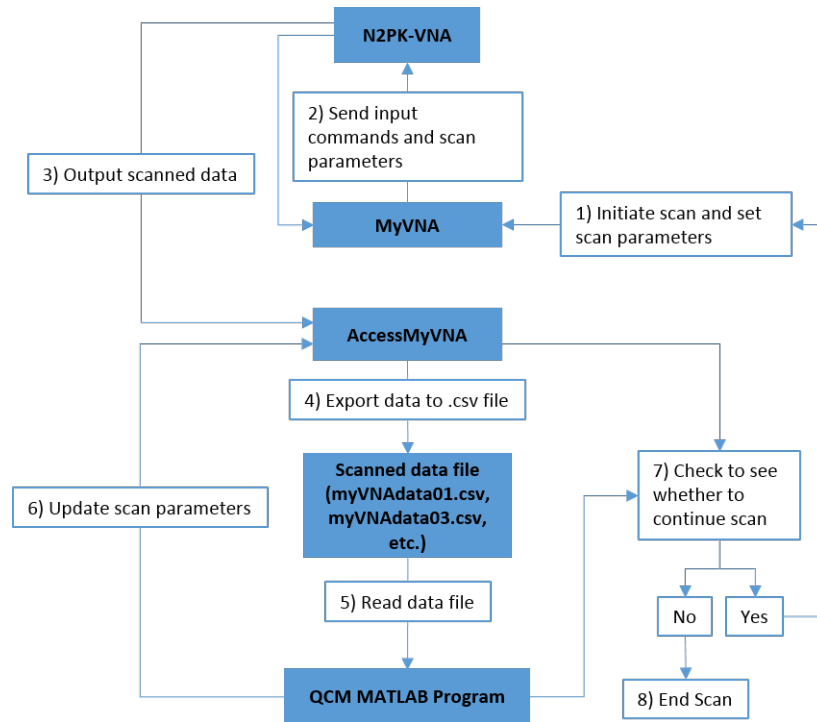


Figure 15: A flow chart of how the three programs communicate with each other.

As shown in Figure 15, the QCM MATLAB Program does not directly communicate with the impedance analyzer, instead the MATLAB program directly manipulates AccessMyVNA to control the impedance analyzer. A more detailed interaction between AccessMyVNA and the QCM MATLAB Program is shown in Figure.

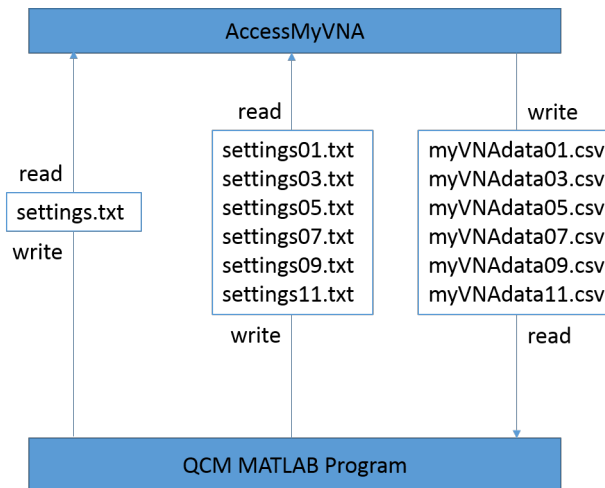


Figure 16: A more detailed flowchart how the QCM MATLAB program communicates with AccessMyVNA.

Part IV.

MATLAB code overview

7. Record scan algorithm explained

Since the program heavily relies on the details on how the measurement process is run, a section has been dedicated in explaining how the MATLAB script was written. When an unforeseen error or bug occurs, this section is a good place in finding the origin of the problem.

8. List of functions

Part V.

Aknowledgements

This program is heavily influenced by the many features in the QTZ.exe written by Diethelm Johannsmann, Institute of Physical Chemistry, TU-Clausthal, Germany. The author of this MATLAB program is very thankful for his contribution to the QCM community and theory. Further, many thanks to the members of the Shull Research Group, Northwestern University, and Professor Kenneth Shull for his guidance and design into the developement of this QCM MATLAB program. Any comments or suggestions are welcome! CHJoshuaYeh@u.northwestern.edu